

Contains Nonbinding Recommendations

Premarket Studies of Implantable Minimally Invasive Glaucoma Surgical (MIGS) Devices

Guidance for Industry and Food and Drug Administration Staff

Document issued on December 15, 2015.

The draft of this document was issued on February 11, 2015.

For questions about this document, contact the Division of Ophthalmic and Ear, Nose, and Throat Devices (DOED) at 301-796-5620.



**U.S. Department of Health and Human Services
Food and Drug Administration
Center for Devices and Radiological Health
Office of Device Evaluation
Division of Ophthalmic and Ear, Nose, and Throat Devices
Intraocular and Corneal Implants Branch**

Preface

Public Comment

You may submit electronic comments and suggestions at any time for Agency consideration to <http://www.regulations.gov> . Submit written comments to the Division of Dockets Management, Food and Drug Administration, 5630 Fishers Lane, Room 1061, (HFA-305), Rockville, MD 20852. Identify all comments with the docket number FDA-2015-D-0288. Comments may not be acted upon by the Agency until the document is next revised or updated.

Additional Copies

Additional copies are available from the Internet. You may also send an e-mail request to CDRH-Guidance@fda.hhs.gov to receive a copy of the guidance. Please use the document number 1400049 to identify the guidance you are requesting.

Table of Contents

I.	Introduction.....	4
II.	Scope.....	5
III.	Definitions	5
IV.	Non-Clinical Testing Recommendations	6
	A. Biocompatibility Testing	7
	B. Physical and Mechanical Testing	9
	C. Sterility and Package Integrity	10
	D. Shelf Life and Shipping Testing.....	11
V.	Clinical Studies	13
	A. Study Design	13
	B. Subject Selection Factors.....	13
	C. Effectiveness Endpoints	15
	D. Safety Outcomes.....	18
VI.	References.....	19

Premarket Studies of Implantable Minimally Invasive Glaucoma Surgical (MIGS) Devices

Guidance for Industry and Food and Drug Administration Staff

This guidance represents the current thinking of the Food and Drug Administration (FDA or Agency) on this topic. It does not establish any rights for any person and is not binding on FDA or the public. You can use an alternative approach if it satisfies the requirements of the applicable statutes and regulations. To discuss an alternative approach, contact the FDA staff or Office responsible for this guidance as listed on the title page.

I. Introduction

This guidance document recommends non-clinical and clinical studies to support a premarket approval (PMA) for implantable minimally-invasive glaucoma surgical (MIGS) devices. Glaucoma is a progressive condition that damages the optic nerve of the eye, is commonly associated with elevated intraocular pressure, and leads to irreversible vision loss. It is the second leading cause of visual disability and blindness in the world. One in 40 adults over 40 years of age has glaucoma severe enough to cause loss of visual function.^{1,2} Current treatments for glaucoma are designed to reduce the intraocular pressure (IOP). Many options are available to lower the IOP including medications, laser treatments, and surgical interventions. Current surgical treatments for glaucoma are aimed at reducing intraocular pressure through the reduction of aqueous inflow or the enhancement of aqueous outflow. While trabeculectomy is the standard surgical intervention for glaucoma, it is often reserved for moderate to severe disease. During the past decade, novel medical devices, called MIGS devices, have emerged. These devices are designed to treat less severe glaucoma by enhancing physiological aqueous outflow with an approach that causes minimal tissue disruption.

This guidance represents the Agency's initial thinking and our recommendations may change as more information becomes available. The Agency strongly encourages manufacturers to engage with CDRH through the Pre-Submission process to obtain more detailed feedback for implantable MIGS devices. For more information on Pre-Submissions, please see "[Requests for Feedback on Medical Device Submissions: The Pre-Submission Program and Meetings with Food and Drug Administration Staff](#)"

Contains Nonbinding Recommendations

(<http://www.fda.gov/downloads/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/UCM311176.pdf>).

For the current edition of the FDA-recognized standard referenced in this document, see the FDA Recognized Consensus Standards Database website at <http://www.accessdata.fda.gov/scripts/cdrh/cfdocs/cfStandards/search.cfm>.

FDA's guidance documents, including this guidance, do not establish legally enforceable responsibilities. Instead, guidances describe the Agency's current thinking on a topic and should be viewed only as recommendations, unless specific regulatory or statutory requirements are cited. The use of the word *should* in Agency guidance means that something is suggested or recommended, but not required.

II. Scope

The recommendations made in this guidance are applicable to implantable MIGS devices, a type of Intraocular Pressure Lowering Implant (associated with product code OGO) used to lower intraocular pressure using an outflow mechanism with either an *ab interno* or *ab externo* approach and associated with little or no scleral dissection and minimal or no conjunctival manipulation. Intraocular Pressure Lowering Implants are Class III devices and are defined as devices intended to reduce intraocular pressure when implanted in eyes which have not failed conventional medical and surgical treatment. This guidance provides additional details on non-clinical and clinical study recommendations not described in the currently FDA-recognized version of ANSI Z80.27: *American National Standard for Ophthalmics – Implantable Glaucoma Devices*.

The recommendations in this guidance document do not apply to implants used to reduce IOP in the anterior chamber of the eye in patients with neovascular glaucoma or with glaucoma when medical or conventional surgical treatments have failed, associated with product code (KYF) and regulated as class II devices under 21 CFR 886.3920, Aqueous Shunt.

III. Definitions

For purposes of this guidance document, the following definitions apply:

Glaucoma: An ophthalmic disease usually characterized by increased intraocular pressure (IOP) resulting in damage to the optic nerve and documented by typical visual field defects.

Humphrey Visual Field (HVF): Automated test method to measure retinal sensitivity. During this test, lights of varying intensities are presented in different parts of the visual field while the subject focuses on one spot. The perception of these lights is charted in decibel (dB) units of light sensitivity.

Hypotony: An intraocular pressure (IOP) less than 6mm Hg.

Contains Nonbinding Recommendations

Intraocular Pressure (IOP): Assessment of pressure in the eye with a tonometer. It is measured in millimeters of mercury (mmHg).

IOP Lowering Implant: A device intended to reduce IOP when implanted in eyes that have not failed conventional medical and surgical treatment.

Hypotony Maculopathy: Abnormality of the macula in the setting of hypotony characterized by optic nerve head swelling, and chorioretinal folds.

Glaucoma Hemifield Test: A particular analysis of the HVF that compares points in the upper field to corresponding points in the lower field and then interprets the results as (a) “outside normal limits” indicating the upper and lower fields are different and may signify glaucoma, (b) borderline, and (c) within normal limits.

Mean Deviation (MD): Average difference between the subject’s overall sensitivity for each point on the visual field when compared with that of the reference population.

Minimally-Invasive Glaucoma Surgical (MIGS) Device: A type of IOP lowering device used to lower IOP using an outflow mechanism with either an *ab interno* or *ab externo* approach, associated with little or no scleral dissection and minimal or no conjunctival manipulation.

Ocular Hypertension: A condition with elevated IOP but no signs of visual field loss or optic nerve damage associated with glaucoma. These subjects may also be called “glaucoma suspects.”

Pattern Deviation (PD) Plot: This measure from the automated visual field provides information about localized defects by adjusting for generalized visual field depression due to non-glaucomatous factors such as media opacity (e.g., cataract or a vitreous hemorrhage).

Washout: Part of a clinical trial when a subject is asked to stop taking all medications. This can occur prior to initiating the investigational treatment as well as before assessing clinical endpoints.

IV. Non-Clinical Testing Recommendations

All non-clinical testing should be performed on the finished sterilized product that is intended to be marketed.

A. Biocompatibility Testing

If the actual device cannot be used in testing (e.g., due to the small area of the device), test samples (e.g., coupons) that are representative of the final device may be employed for biocompatibility testing.

1. Recommended Biocompatibility Tests

The following tests should be performed as recommended by [Bluebook Memorandum G95-1 Use of International Standard ISO-10993, “Biological Evaluation of Medical Devices Part 1: Evaluation and Testing.”](#) (<http://www.fda.gov/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/ucm080735.htm>).

- a. Cytotoxicity
- b. Sensitization
- c. Ocular irritation
- d. Systemic toxicity (acute toxicity)
- e. Sub-chronic toxicity (subacute toxicity)
- f. Genotoxicity
- g. Carcinogenicity
- h. Pyrogens Testing. If the device contacts blood then material-mediated pyrogenicity testing is also recommended.

In addition, ocular implantation testing should be conducted as outlined in Annex B of the currently FDA-recognized version of ANSI Z80.27: *American National Standard for Ophthalmics – Implantable Glaucoma Devices*. There might be cases (e.g., inflammation) in which the 6-month implantation study recommended in ANSI Z80.27 is not sufficient and longer implantation studies may be needed. For recommended implant test periods for various animal models, please refer to FDA-recognized sections of ISO 10993-6: *Biological evaluation of medical devices — Part 6: Tests for local effects after implantation*.

2. Recommended Physico-Chemical Tests

- a. Test of Extractables and Hydrolytic Stability: Testing should be conducted as outlined in Annex C of the currently FDA-recognized version of ISO 11979-5: *Ophthalmic Implants – Intraocular Lenses – Part 5: Biocompatibility*.
- b. Test of Extractables by Exhaustive Extraction (Annex C of ISO 11979-5)
- c. Leachables (Annex D of ISO 11979-5)
- d. Insoluble Inorganics (ISO 11979-5)

Contains Nonbinding Recommendations

3. Bioabsorbable Materials

This testing should be performed if the material is in situ polymerizing and bioabsorbable.

Toxicity should be assessed for the finished product as well as at various time points over the course of polymerization and/or degradation to ensure that starting, intermediate and final degradation products are evaluated. Assessments should continue until the polymer is no longer present in the tissue, or until the biological tissue response is demonstrated to be stable. For additional guidance, please refer to the currently FDA-recognized sections of ISO 10993-6 and the currently FDA-recognized version of ANSI Z80.27.

4. Biological Response from Device Mechanical Failure

For devices incorporating a coating or multiple material components, it is possible that mechanical failure could alter the biological response to the device. For devices with the potential for biological hazard due to mechanical failure, the biocompatibility testing should include testing to address this concern. We encourage manufacturers to propose a minimum limit for device failure as determined with scientific justification. This should be based on considerations such as the device material and design, risk analysis, anatomical location, and insertion method.

5. Sample Preparation

For biocompatibility testing using extracts of samples, the extraction should be conducted using both polar (water, physiological saline) and non-polar (sesame oil, cotton oil) extraction vehicles under conditions as described in the currently FDA-recognized version of ISO 10993-12: *Biological evaluation of medical devices -- Part 12: Sample preparation and reference materials*. For permanently implanted devices, extraction at 37°C for 72 hours may not be sufficient to obtain an extract that represents the chemicals that may leach out over the use life of the device. However, in some cases, temperatures over 37°C may result in degradants and toxicities that are not representative of the device. Therefore, a justification for the selected extraction conditions should be provided.

Extraction should be performed based on surface area of the device. If the area cannot be determined than a mass/volume should be used. The test extract should not be processed (e.g., filtered or centrifuged) and should be used immediately after preparation.

Extraction in tissue culture media supplemented with serum is acceptable for cytotoxicity testing and should be performed according to the currently FDA-recognized version of ISO 10993-5: *Biological evaluation of medical devices -- Part 5: Tests for in vitro cytotoxicity*.

Contains Nonbinding Recommendations

A scientifically-based rationale for omission of any recommended test should be included with the submission. We recommend that sponsors who do not intend to conduct biocompatibility testing submit a pre-submission to obtain feedback from the Division of Ophthalmic and Ear, Nose, and Throat Devices on their rationale. For more information on Pre-Submissions, please see FDA's guidance "[Medical Devices: The Pre-Submission Program and Meetings with FDA Staff](http://www.fda.gov/downloads/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/UCM311176.pdf)" (<http://www.fda.gov/downloads/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/UCM311176.pdf>).

B. Physical and Mechanical Testing

Device properties should be determined at *in situ* conditions with the temperature tolerance of ± 2 °C. The precise composition of the solution used should be reported in all cases. FDA recommends that testing be conducted as outlined in Physical and Mechanical Testing of Section 5 of ANSI Z80.27 with the following additions and exceptions.

1. Validation of Dimensional Tolerances

(Section 5.4 of ANSI Z80.27) Dimensions for which tolerances are given should be specified in the manufacturer's design documentation. The sponsor should validate that their production meets their tolerances to appropriate statistical levels.

2. Surface and Edge Quality

(Sections 5.2 and 5.3 of ANSI Z80.27) The device should be essentially free from surface defects and all edges should appear smooth when viewed at 10x magnification with a stereo microscope using optimal lighting conditions. Any questionable or critical areas should be viewed at higher magnification.

3. Structural Integrity

(Section 5.7 of ANSI Z80.27) The manufacturer should provide evidence that the device can withstand surgical manipulations without failure. An appropriate test method and specification should be established by the manufacturer to ensure that the device does not fail at typical deformations.

4. Insertion Testing

The purpose of this testing is to evaluate the integrity of the delivery system and of the delivered device, if the MIGS device is designed to be delivered from an injector system. The injector system should be evaluated following the instructions supplied by the manufacturer and using recommended lubricants and

Contains Nonbinding Recommendations

instrumentation. There should be no change in the physical properties of the MIGS device and no damage to the injector system as a result of the delivery. The results should be reported and are acceptable if the physical properties of the MIGS device remain within manufacturing specifications of the product.

5. Coated Devices

MIGS devices with surface coatings should conduct testing per Section 9.2 of ANSI Z80.27.

6. Metallic Devices

MIGS devices manufactured with metallic materials should be evaluated for Magnetic Resonance Imaging (MRI) safety according to “[FDA Guidance for Industry and FDA Staff: Establishing Safety and Compatibility of Passive Implants in the Magnetic Resonance \(MR\) Environment](http://www.fda.gov/downloads/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/UCM107708.pdf)” (<http://www.fda.gov/downloads/MedicalDevices/DeviceRegulationandGuidance/GuidanceDocuments/UCM107708.pdf>) and for corrosion resistance according to the currently FDA-recognized version of ASTM F2129: *Standard Test Method for Conducting Cyclic Potentiodynamic Polarization Measurements to Determine the Corrosion Susceptibility of Small Implant Devices*.

C. Sterility and Package Integrity

1. Sterilization Method

The sterilization method should be validated according to one of the following standards:

- a. For moist heat (steam), use the currently FDA-recognized version of ANSI/AAMI/ISO 17665-1: *Sterilization of Health Care Products – Moist Heat – Part 1: Requirements for the Development, Validation, and Routine Control of a Sterilization Process for Medical Devices*.
- b. For ethylene oxide, use the currently FDA-recognized version of ISO 11135: *Sterilization of Health Care Products – Ethylene Oxide – Requirements for the Development, Validation, and Routine Control of a Sterilization Process for Medical Devices*.
- c. For gamma radiation, use the currently FDA-recognized version of ANSI/AAMI/ISO 11137-1: *Sterilization of Health Care Products – Radiation – Part 1: Requirements for the Development, Validation, and Routine Control of a Sterilization Process for Medical Devices*.

Contains Nonbinding Recommendations

2. Ethylene Oxide Sterilant Residues

If the MIGS device is sterilized via ethylene oxide, then the maximum level of ethylene oxide residuals that remain on the device should be quantified and assessed according to the currently FDA-recognized version of ISO 10993-7: *Biological evaluation of medical devices – Part 7: Ethylene oxide sterilization residuals*. An exhaustive solvent or head space extraction method should be chosen and the amount of residue should conform to those for intraocular lenses. If the extraction is not exhaustive, release criteria should be lowered in proportion to the relative efficiency of the method.

The residue of ethylene chlorohydrin should not exceed a release of more than 2.0 µg per device per day and not exceed 5.0 µg in total per device.

3. Bacterial Endotoxins

The recommended endotoxin limit for MIGS devices is ≤ 0.2 EU/device as stated in the FDA's guidance "[Endotoxin Testing Recommendations for Single-Use Intraocular Ophthalmic Devices](http://www.fda.gov/downloads/medicaldevices/deviceregulationandguidance/guidancedocuments/ucm393376.pdf)" (<http://www.fda.gov/downloads/medicaldevices/deviceregulationandguidance/guidancedocuments/ucm393376.pdf>). This limit applies to the segment of the device placed in the anterior chamber and the segment(s) contacting the aqueous humor even though the main portion of the device may reside outside the eye. For devices that have a segment that contacts the aqueous humor and the vitreous or posterior segment, please contact the Division.

4. Package Integrity Testing

Package integrity testing should be performed regardless of the sterilization method and may consist of a validated whole package physical integrity test in combination with a validated seal integrity test. Examples of whole package physical integrity testing can be found in FDA's guidance "[Container and Closure System Integrity Testing in Lieu of Sterility Testing as a Component of the Stability Protocol for Sterile Products](http://www.fda.gov/RegulatoryInformation/Guidances/ucm146074.htm)" (<http://www.fda.gov/RegulatoryInformation/Guidances/ucm146074.htm>) or the currently FDA-recognized version of ANSI/AAMI/ISO 11607-1: *Packaging for Terminally Sterilized Medical Devices – Part 1: Requirements for Materials, Sterile Barrier Systems and Packaging Systems*.

D. Shelf Life and Shipping Testing

1. Development of Shelf Life Protocol

The protocol for the shelf life study should be developed prior to initiation of the study. If, during the course of the study, a parameter no longer conforms to the

Contains Nonbinding Recommendations

manufacturing specifications at two or more time intervals, the maximum shelf-life of the MIGS device under study has been reached at the last conforming measurement point. If a manufacturer wishes to maintain the possibility to re-sterilize finished device lots, the finished device lot(s) used in the stability study should undergo the maximum number of sterilization cycles allowed under the manufacturer's procedures. References to suggested test methods can be found in the currently FDA-recognized version of ISO 11979-6: *Ophthalmic Implants – Intraocular Lenses – Part 6: Shelf-life and Transport Stability*.

2. Real-time Shelf-Life Study

FDA recommends conducting the following stability and integrity studies:

- a. Product Stability Studies**
 - (1) Dimensions
 - (2) Surface and Edge Quality
 - (3) Structural Integrity
 - (4) Pressure/Flow Characterization
 - (5) Insertion Testing
 - (6) Coating Stability, if applicable

- b. Package Integrity Studies**
 - (1) Whole Package Physical Integrity
 - (2) Seal/Closure Integrity

3. Accelerated Shelf-Life Studies

These studies are the same as those performed for real-time shelf life studies with the exception of the conditions in which they are performed. It is important that devices to be measured be allowed to equilibrate to the same conditions as at the initial measurements before being tested. The corresponding real-time shelf-life is calculated by multiplying the studied time period by $2^{(T_a - T_o)/10}$, where T_a is the accelerated temperature and T_o is the typical storage temperature (usually room temperature). The maximum acceptable storage temperature is 45°C. While an initial shelf-life can be established with accelerated testing, a confirmatory real-time shelf-life study should be performed.

4. Transport Stability

The complete, filled device packages (in their normal transport package) should be able to withstand extremes of the temperature and humidity (as expected in shipping), vibration and being dropped. Both the packaging and the product should be inspected following completion of the pre-test conditioning. The device should be considered to have satisfactorily passed the test if the device is free from physical damage when visually inspected under magnification. The packaging should also continue to provide functional protection to the device.

Contains Nonbinding Recommendations

FDA recommends that the following tests be performed, at a minimum:

- a. Legibility of Labeling (empty packages can be used);
- b. Surface and Edge Quality (sealed packages should be used);
- c. Seal/Closure Integrity (empty packages can be used);
- d. Whole Package Physical Integrity (empty packages can be used).

V. Clinical Studies

A. Study Design

1. Patient Follow-up

It is strongly recommended that all subjects be followed for a minimum of 12 months prior to submission of any premarket application, as discussed at the FDA/AGS Workshop on Supporting Innovation for Safe and Effective Minimally Invasive Glaucoma Surgery, February 26, 2014. For additional information, refer to the workshop materials and transcript available on FDA's website at <http://www.fda.gov/MedicalDevices/NewsEvents/WorkshopsConferences/ucm382508.htm>. For follow-up of less than 24 months, you should provide justification based upon the benefit-risk analysis. For further information on the principal factors FDA considers when making benefit-risk determinations during the premarket review process, please see "[Factors to Consider When Making Benefit-Risk Determinations in Medical Device Premarket Approvals and De Novo Classifications](http://www.fda.gov/medicaldevices/deviceregulationandguidance/guidancedocuments/ucm267829.htm)" (<http://www.fda.gov/medicaldevices/deviceregulationandguidance/guidancedocuments/ucm267829.htm>). If the benefit-risk analysis raises concerns beyond 24 months after implantation, longer follow-up may be appropriate. The investigational plan should include the possibility that long-term follow-up (e.g., up to five years) may be necessary. It is recommended that informed consent for up to five years of follow-up is obtained.

2. Washout Period

All eligible subjects should undergo a washout prior to surgery and at the key timepoints of interest for the study.

B. Subject Selection Factors

Subjects included in clinical trials for MIGS devices should have evidence of early or moderate open angle (confirmed by gonioscopy) glaucoma.* For the purpose of this

* Based on discussion at the FDA/ AGS Workshop on Supporting Innovation for Safe and Effective Minimally Invasive Glaucoma Surgery, February 26, 2014, we do not believe that pre-perimetric glaucoma (i.e., optical coherence tomography (OCT) changes and optic nerve changes without any field abnormalities) should be included

Contains Nonbinding Recommendations

guidance, glaucoma is defined as characteristic visual field and optic nerve changes as outlined below.

1. Inclusion Criteria

a. Humphrey Visual Field (HVF)

The HVF should be reliable, which is defined as fixation losses, false positives, and false negatives all less than 33%.³ The following characteristics should also be noted on the HVF:

- (1) Visual field defects consistent with glaucomatous optic nerve damage;⁴ and
- (2) Mean deviation not worse than -12 dB; and at least one of the following two findings:
 - (a) A cluster of 3 or more points in an expected location of the visual field depressed below the 5% level, at least 1 of which is depressed below the 1% level on the pattern deviation (PD) plot
 - (b) Glaucoma hemi-field test “outside normal limits.”

While these recommendations are based on HVF, equivalent perimetric criteria using a different device could be considered for trial enrollment.

b. Glaucomatous Optic Nerve Damage

Glaucomatous optic nerve damage as evidenced by any of the following optic disc or retinal nerve fiber layer structural abnormalities:

- (1) Diffuse thinning, focal narrowing, or notching of the optic disc rim, especially at the inferior or superior poles with or without disc hemorrhage;
- (2) Localized abnormalities of the peripapillary retinal nerve fiber layer, especially at the inferior or superior poles; or
- (3) Optic disc neural rim asymmetry of the two eyes consistent with loss of neural tissue

2. Exclusion Criteria

Subjects that should be excluded from clinical trials for MIGS devices include but are not limited to the following:

in these interventional studies because there are differing opinions amongst experts about whether this condition warrants surgical treatment. Because this guidance focuses on an emerging topic, it represents the Agency’s initial thinking and our recommendations may change as more information becomes available. The current recommendations are designed to provide a conservative approach to protection of human subjects.

Contains Nonbinding Recommendations

- a.** Subjects that cannot perform HVF or equivalent perimetric testing in either eye
- b.** Subjects who cannot undergo a medication “washout” in the study eye or who are at high risk for adverse outcomes, including subjects:
 - (1) Requiring oral medications (e.g., acetazolamide) for IOP control
 - (2) With severe glaucoma in either eye, defined as mean deviation (MD) of -12.00 or worse and at least one of the following:
 - (a) On PD plot, greater than or equal to 75% of points depressed below the 5% level and greater than or equal to 50% of points depressed below 1% level; or
 - (b) At least 50% of points (i.e., 2 or more) within central 5 degrees with sensitivity of < 0dB on the dB plot; or
 - (c) Points within the central 5 degrees of fixation with sensitivity < 15dB in both hemifields on the dB plot.
 - (3) With fixation-threatening glaucoma in either eye noted at the qualifying visit: Subjects with visual field defects threatening fixation defined as any (1 or more) point(s) within the central 5° depressed below the 5th percentile on PD plot unless this/these points are >25 dB on Threshold Values (decibel) plot.
 - (4) With best corrected visual acuity of 20/200 or worse in the study eye due to glaucoma.
- c.** Subjects with ocular hypertension
- d.** Subjects at high risk for adverse outcomes due to placing a device in the angle

For details of other subject inclusion and exclusion characteristics (e.g., minimal endothelial cell density), please refer to the non-refractory section of the currently FDA recognized version of ANSI Z80.27: *American National Standard for Ophthalmics – Implantable Glaucoma Devices*.

C. Effectiveness Endpoints

1. Washout

All subjects should undergo a washout period of all IOP-lowering medications prior to surgery to establish a baseline IOP. In addition, if IOP-lowering medications are re-instituted postoperatively, all subjects should undergo a

Contains Nonbinding Recommendations

washout period prior to the time point(s) for data collection used in the effectiveness analyses.

2. Primary effectiveness

The recommended primary effectiveness endpoint is the percentage of subjects with reduction of at least 20% (i.e., $\geq 20\%$) in mean diurnal IOP from baseline.⁵⁻⁹ The proposed hypothesis test for the primary effectiveness endpoint should be described in the statistical analysis plan.

3. Secondary effectiveness

The recommended secondary effectiveness endpoint is the mean diurnal IOP change from baseline.

4. Recommended Analyses

In addition to the analyses described in ANSI Z80.27: *American National Standard for Ophthalmics – Implantable Glaucoma Devices*, FDA recommends the following:

a. Scatter Plots of Postoperative IOP as a function of Preoperative IOP

Scatter plots should be provided for the mean diurnal IOP measurements. Reference lines should be added to the scatter plots at various postoperative IOP levels of interest (e.g., 18 mmHg), and the percentage of subjects falling on either side of each of the reference lines should be presented (e.g., the percent of subjects achieving postoperative IOP of less than 18 mmHg and the percent with postoperative IOP of greater than or equal to 18 mmHg). Appropriate reference lines should be selected based upon such considerations as the specific enrollment criteria for the study, such as the IOP and severity of disease criteria.

b. Percent Reduction in Mean Diurnal IOP

The number and percent (e.g., n/N & %) of subjects achieving percent reduction (or increase) in mean diurnal IOP at each annual visit from baseline across the entire distribution should be presented. This analysis should be presented with and without further stratification by baseline mean diurnal IOP.

c. Changes in the Mean, Range, and Maximum of the Diurnal IOP Measurements, and Box-plots of Mean, Range, and Maximum of Diurnal IOP Measurements

Contains Nonbinding Recommendations


In addition to analyses of the changes in the mean diurnal IOP measurements from baseline described in ANSI Z80.27, , analyses of the changes in the range of the diurnal IOP measurements and maximum of the diurnal IOP measurements from baseline should be conducted. Descriptive statistics should be performed as described in ANSI Z80.27 with additional stratification by baseline mean diurnal IOP. In addition to the box plots of mean diurnal IOP measurements outlined in the standard, box plots of the range of diurnal IOP measurements and of the maximum diurnal IOP measurements should be presented. An example of a box-plot representation of IOP data can be found in the World Glaucoma Association's (WGA) Guidelines on Design and Reporting of Glaucoma Surgical Trials, Consensus on definitions of success, Section II, General Data Presentation Requirements.¹⁰

d. Fluctuation of IOP Measurements Over Time

For each subject, we recommend plotting the diurnal IOP measurements (y-axis) versus time of measurements (x-axis) for baseline and each of the postoperative diurnal IOP visits on the same graph using a different symbol for each visit (See examples in Figures 1 and 2).¹⁰

e. Change in Number of Medications

At each postoperative visit, a graphical representation of the number pre-operative (before washout, when applicable) on the x-axis versus post-operative IOP-lowering medications (counting combination drops as separate medications) on the y-axis should be made. An example of such a graphical representation is presented below in Figure 3. The size of each bubble represents the number of subjects.

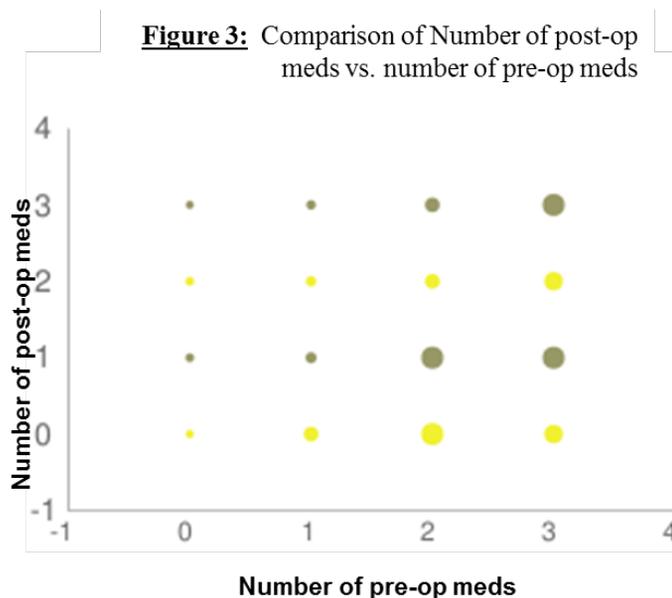

Contains Nonbinding Recommendations

f. Assessment of Balance in Baseline Variables

In addition to the analyses recommended in ANSI Z80.27 and above, we recommend checking for imbalances in baseline variables that may affect the outcome (e.g., baseline IOP, age, race, gender, number of medications at screening) among the arms of the trial.

D. Safety Outcomes

1. Adverse Events and Device Malfunctions

The adverse events and device malfunctions for MIGS devices are listed in ANSI Z80.27. The definition of each adverse event should specify the grade or severity, the degree of involvement of the anatomical structure, the timing, and the duration of the event, as applicable, in order to distinguish findings that should be reported as “adverse events” from those observations that should be routinely recorded. Case report forms should include a forced-choice method of recording listed adverse events as well as a method of recording other adverse events not listed.

- a.** Hypotony should be classified as an early (i.e., at 2 weeks or less following surgery) or late (i.e., more than 2 weeks after surgery) adverse event if it occurs with at least one of the following conditions:
 - (1) Flat anterior chamber requiring anterior chamber reformation
 - (2) Corneal folds
 - (3) Choroidal effusions requiring or undergoing surgical drainage
 - (4) Suprachoroidal hemorrhage
 - (5) Maculopathy
 - (6) Corneal astigmatism

- b.** Chronic anterior uveitis should be defined as inflammation of grade 1+ or worse persisting for more than 3 months post-operatively or that recurs less than three months after discontinuing treatment.¹¹ Two line or greater decrease in best corrected visual acuity (BCVA)

2. Device Removal

Reasons for device removal should be recorded in the case report forms.

VI. References

1. Resnikoff S, Pascolini D, Etya'ale D, et al. Global data on visual impairment in the year 2002. *Bulletin of the World Health Organization*. 2004;82(11):844-851.
2. Quigley HA. Glaucoma. *Lancet*. 2011;377(9774):1367-1377.
3. Birt CM, Shin DH, Samudrala V, Hughes BA, Kim C, Lee D. Analysis of reliability indices from Humphrey visual field tests in an urban glaucoma population. *Ophthalmology*. 1997;104(7):1126-1130.
4. Hodapp E, Parrish RK, 2nd, Anderson DW. *Clinical decisions in glaucoma*. St. Louis, MO: The CV Mosby Co; 1993.
5. Craven ER, Katz LJ, Wells JM, Giamporcaro JE, iStent Study G. Cataract surgery with trabecular micro-bypass stent implantation in patients with mild-to-moderate open-angle glaucoma and cataract: two-year follow-up. *J Cataract Refract Surg*. 2012;38(8):1339-1345.
6. Gedde SJ, Schiffman JC, Feuer WJ, et al. The tube versus trabeculectomy study: design and baseline characteristics of study patients. *Am J Ophthalmol*. 2005;140(2):275-287.
7. Barton K, Gedde SJ, Budenz DL, Feuer WJ, Schiffman J, Ahmed Baerveldt Comparison Study G. The Ahmed Baerveldt Comparison Study methodology, baseline patient characteristics, and intraoperative complications. *Ophthalmology*. 2011;118(3):435-442.
8. Christakis PG, Tsai JC, Zurakowski D, Kalenak JW, Cantor LB, Ahmed, II. The Ahmed Versus Baerveldt study: design, baseline patient characteristics, and intraoperative complications. *Ophthalmology*. 2011;118(11):2172-2179.
9. Samuelson TW, Katz LJ, Wells JM, Duh YJ, Giamporcaro JE. Randomized evaluation of the trabecular micro-bypass stent with phacoemulsification in patients with glaucoma and cataract. *Ophthalmology*. 2011;118(3):459-467.
10. Heuer DK, Barton K, Grehn F, Shaarawy TM, Sherwood MB. Consensus on definitions of success. In: Shaarawy TM, Sherwood MB, Grehn F, eds. *World Glaucoma Association Guidelines on Design and Reporting of Glaucoma Surgical Trials*. Amsterdam: Kugler Publications; 2009.
11. Jabs DA, Nussenblatt RB, Rosenbaum JT, Standardization of Uveitis Nomenclature Working G. Standardization of uveitis nomenclature for reporting clinical data. Results of the First International Workshop. *Am J Ophthalmol*. Sep 2005;140(3):509-516.